
BT-CSIS

Tube Feeding Training Model

USER'S GUIDE

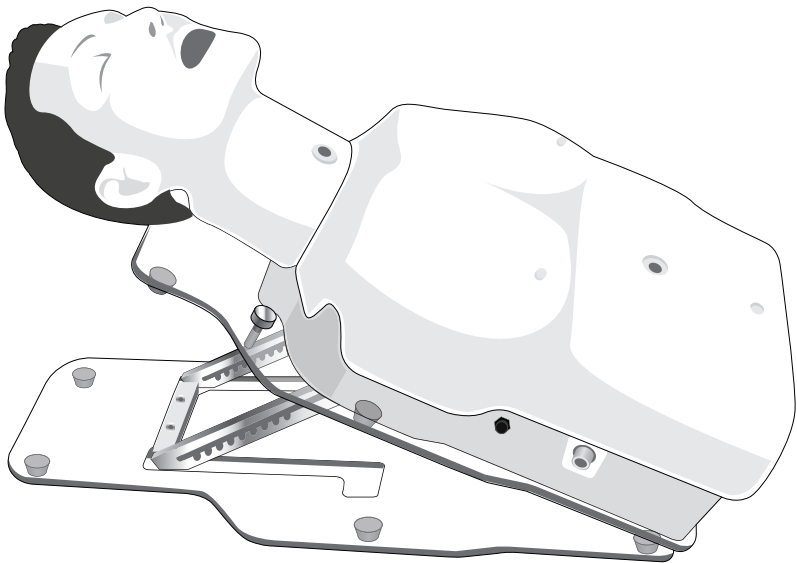■ Product Information

Product Name	Tube Feeding Training Model
Model Name	BT-CSIS
Manufacturer	BT Inc.
Country of Origin	South Korea

■ Product Specifications

Size	310 x 700 x 238 mm
Weight	9kg

■ Customer Service

Web-site	www.btinc.co.kr
Address	Headquarter : A-313 Samsong Techno-valley, 140, Tongil-ro, Deogyang-gu, Goyang-si, Gyeonggi-do, 10594, South Korea
Contacts	Sales : Busan Building 5F, 27, Jeonpo-daero 92beon-gil, Nam-gu, Busan 48404, South Korea
E-mail	Headquarter : +82-2-2038-2903 Sales : +82-51-903-1144 bt@btinc.co.kr



I Cautions

Please read it carefully to ensure correct usage and user's safety through understanding of the product.

- ① Use only components provided or approved by BT Inc.
- ② Do not use the product for other purposes described in the manual.
- ③ If you disassemble or modify the product without any permission by BT Inc, you cannot get the warranty. Also, making imitated or altered products for commercial purposes are prohibited by law.
- ④ Keep the product away from water and moisture.
- ⑤ Do not keep the product near hot or cold environments as possible. Deterioration of product or internal damage may occur by the extreme weather conditions
- ⑥ Never clean the product with chemicals such as solvent, benzene or thinner as these may cause fire ignition or the surface deterioration.
- ⑦ BT Inc. can change or adjust the product specification without any notice in order to improve the product.

Components

Adult torso model

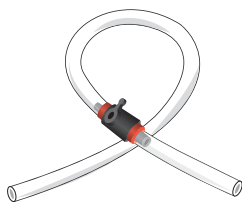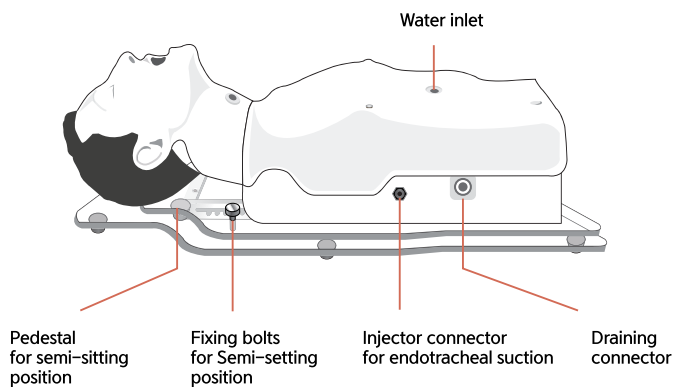Draining tube

Funnel



Injector for endotracheal suction



Inner stopper for stomach

* Etc. – Storage bag, User manual



Installation

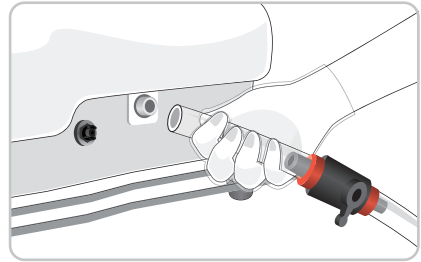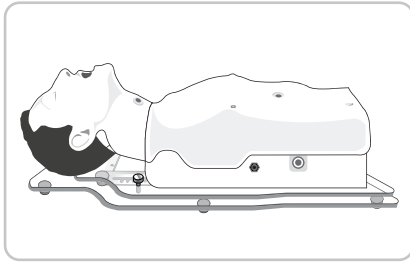Put the manikin in a supine position, lock the valve of the draining tube and insert it into the draining connector of manikin.

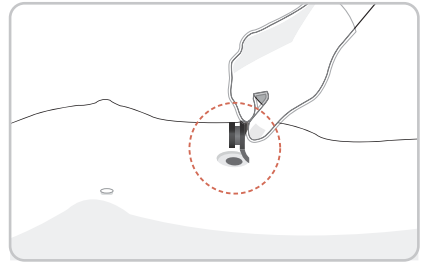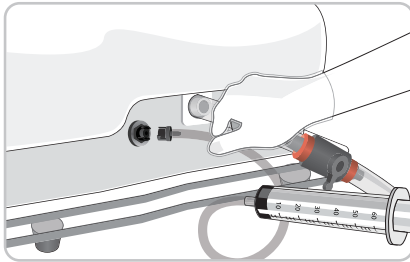Insert the injector for endotracheal suction to the connector and open the valve of water inlet on the abdomen.

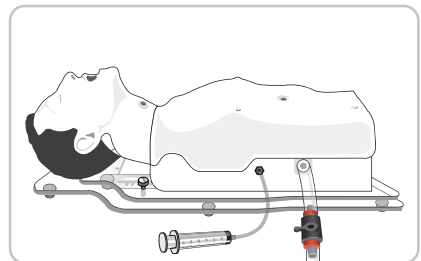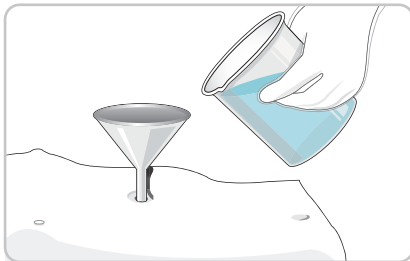
Remove the inner stopper inside the stomach, fill about 300cc ~ 350cc of water by using a funnel. After that, insert the inner stopper in the stomach part and close the outer cover.


※ When install the tube feeding training model,
be sure to lock the valve of draining tube before injecting water.

Installation


<Changing position>

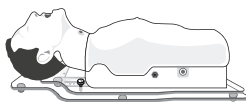Supine position


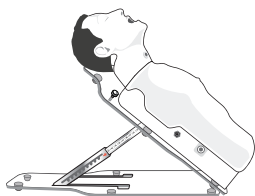
Semi-sitting position (21°~45°)

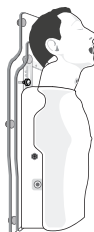
Sitting position

<How to change the Semi-sitting position>

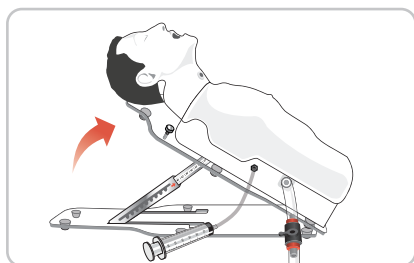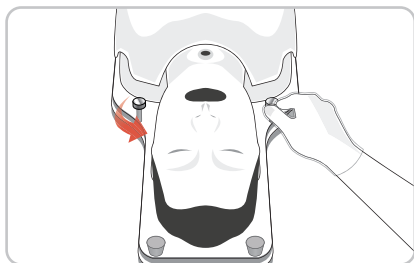
Loosen and remove bolts at the top, adjust the pedestal to the angle.

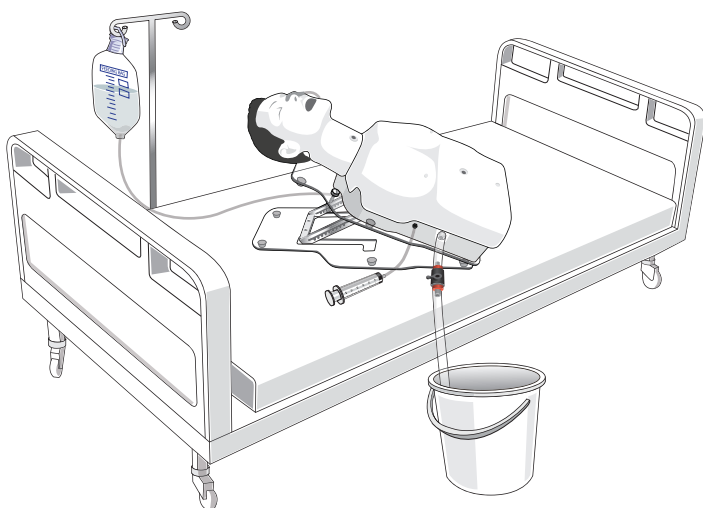
Functions

<Tube feeding>

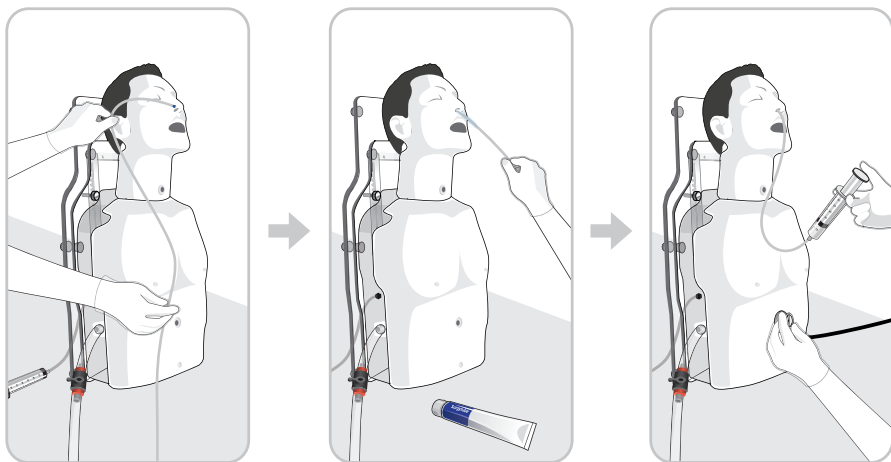When practicing tube feeding, set the manikin to sitting or semi-sitting position and practice according to the way of inserting the catheter.

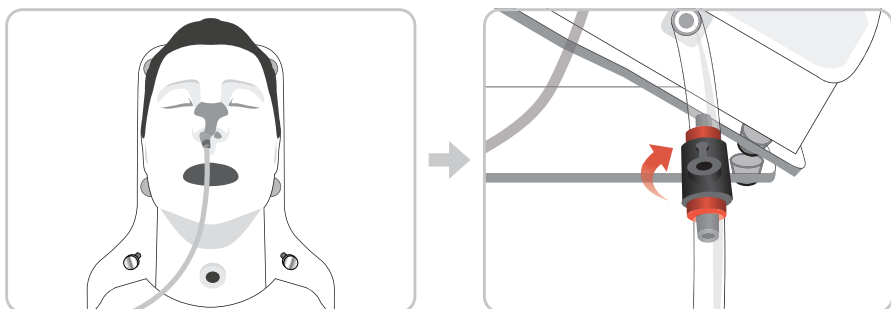
When injecting nutrition after inserting a catheter, inject nutrition after open the valve of darining tube.

※ When practicing tube feeding, be sure to open the valve of draining tube.

Functions


<Tracheostomy care / Endotracheal suction>

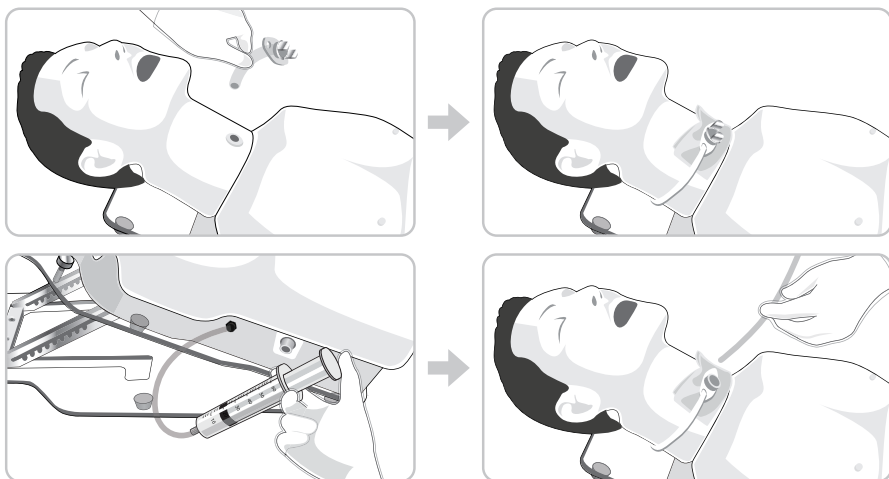Practicing tracheostomy care is possible.


Please do the training after injecting artificial saliva of 20~30cc with a injector for endotracheal suction.

<Gastrostomy tube / PEG care>

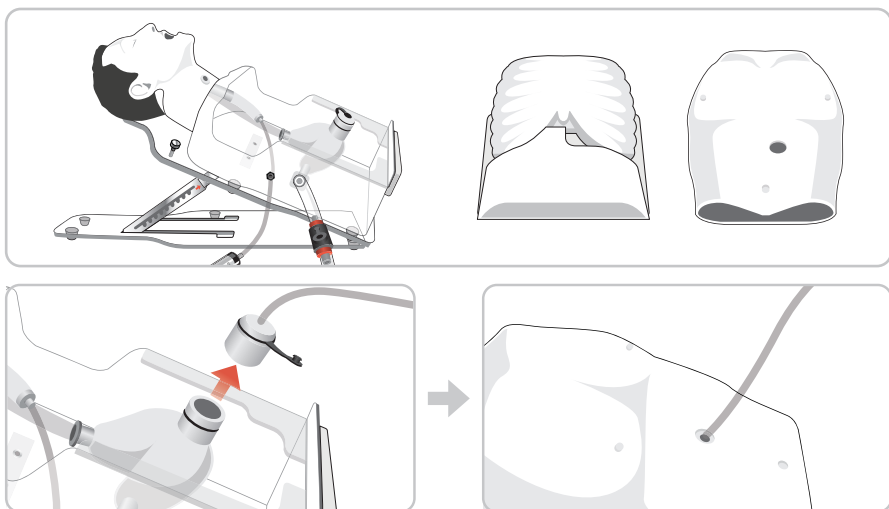
When practicing gastrostomy tube care or PEG care, remove the chest skin and rib part of manikin, insert a gastrostomy tube by opening the cap of stomach, seal the skin and proceed the training.


Maintenance

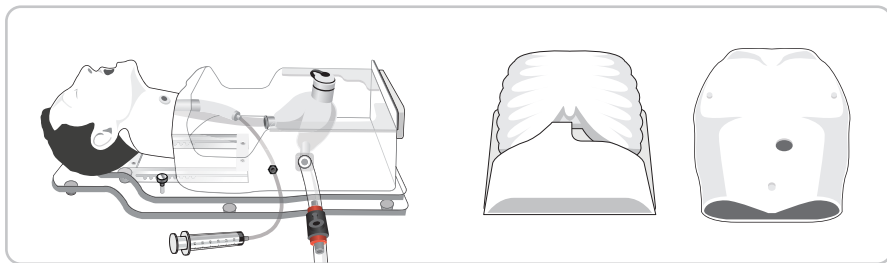After training, put the manikin in a lateral position, drain the remaining nutrition fluid and remove the chest skin and rib part of manikin.

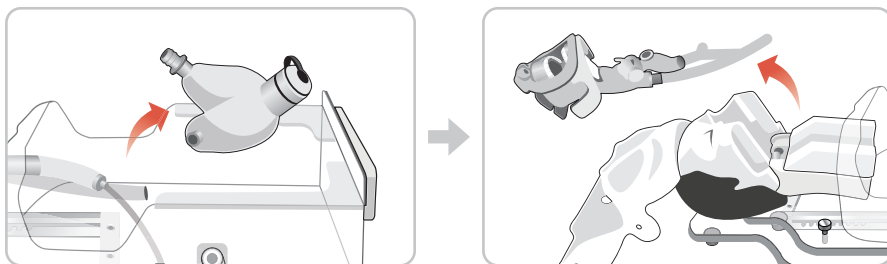Separate the stomach, remove the face skin, and then separate the trachea.

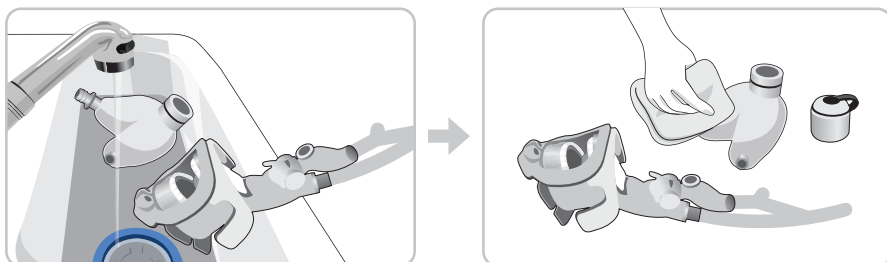Wash models in running water, wipe the surface with a dry cloth and dry it. After drying, assemble and store in reverse order.

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Copyright(c) BT Inc. All rights reserved.