

# INSTRUCTION MANUAL

## Stump Bandaging Simulators

LF01063 & LF01064



# About the Simulator



Featuring the realism of appearance always found with Nasco Healthcare Life/form® products, these simulators consist of two separate parts to make the learning and demonstrating of stump bandaging much more effective and convenient for teacher and student. These simulators eliminate the need for live patients and gives the students an opportunity to practice at their own pace.

# List of Components

The upper torso includes two arms, one amputated above the elbow, and the other above the wrist. Both arms are slightly extended, as a patient would have their arm while being bandaged. The compressibility of the material very closely duplicates that of a live patient and makes the experience very lifelike. The lower half is designed in the shape of a prone patient with legs slightly abducted. One leg is amputated below the knee, and the other above the ankle. The stumps are formed to add to the realism.

## Procedures That Can Be Performed On These Simulators

Most standard bandaging procedures can be practiced or demonstrated. The lower torso is long enough to allow carrying the bandage around the body over the uninvolvement hip at the level of the iliac crest. (Note: There is no embedded palpable skeletal material.) Both units can be used to demonstrate the attachment of prosthetic devices.

## Care Of Simulator

Normal soil accumulated on the surface of the simulator can be removed with mild soap and lukewarm water. Use Nasco Cleaner (LF09919) to remove stubborn stains from the simulator. Simply spray soiled area and wipe clean with a soft cloth or paper towels. The most kinds of adhesive bandages can easily be removed. Use only Glue-Off (LF01023) for removing residual adhesives. Apply a thin coat. Allow to dry, then rub off in a circular motion.

## Cautions


Solvents or corrosive materials will damage the simulator. Never place simulator on any kind of printed paper or plastic. These materials will transfer indelible stains. Ballpoint pens will also make indelible stains.

## Supplies/Replacement Parts For Stump Bandaging Simulators

LF01023 Glue-Off  
LF09919 Nasco Cleaner

# Other Compatible Nasco Healthcare Simulators

- LF00929U** Surgical Bandaging
- LF01028U** Suture Practice Arm
- LF01034U** Suture Practice Leg

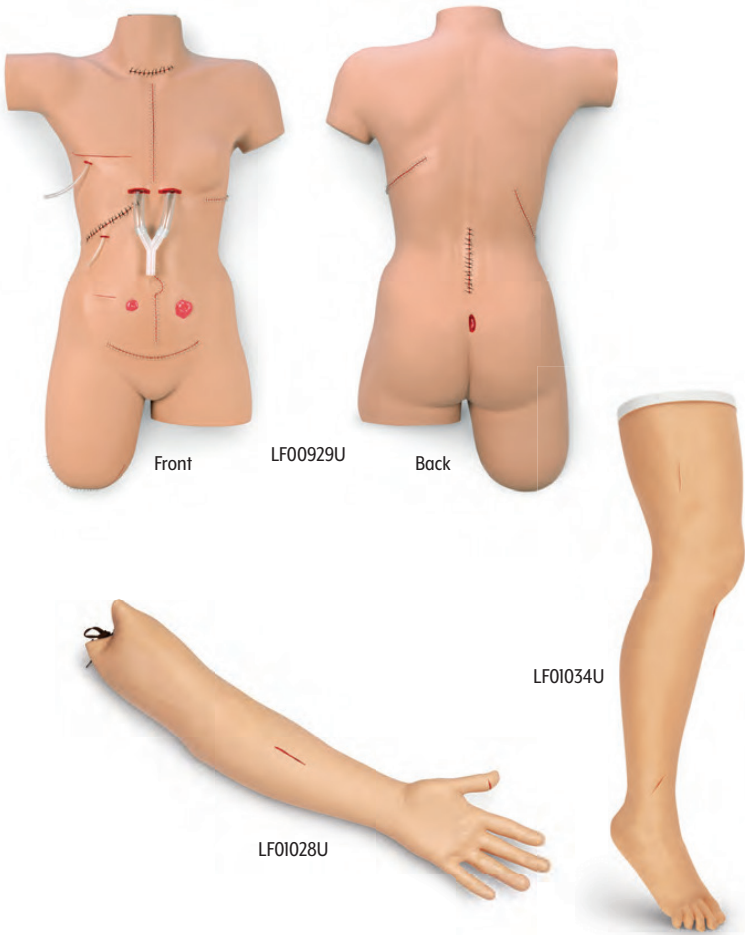

**Nasco Healthcare**  
16 Simulaids Drive  
Saugerties, NY 12477  
1-833-NASCOHC (627-2642)  
info@nascohealthcare.com  
www.nascohealthcare.com